


# ***Gamma-2 Abyss Deep-Sea Transmissometer***

## **USER'S MANUAL**

Revision E

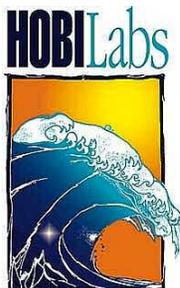

---

**Hydro-Optics, Biology & Instrumentation  
Laboratories**

*LIGHTING THE WAY IN AQUATIC SCIENCE*

[www.hobilabs.com](http://www.hobilabs.com)

[support@hobilabs.com](mailto:support@hobilabs.com)

Revisions:

E, April, 2011: Add /S switch to BAUD command (section 5.2.1); add argument options to DIR command (section 5.2.6); add /N switch to START command (section 5.2.12); add VSLEEP command (section 5.2.18); add section 7.

D, March, 2011: Change “Abyss-2” to “Gamma-2 Abyss.”

C, November 2010: Correct calibration formulas (section 6.4.2)

B, May 2010: Extensive revision and additions.

A, May 2010: Initial Release

<b>1</b>	<b>INTRODUCTION .....</b>	<b>4</b>
1.1	OPTO-ELECTRONIC DESIGN .....	4
1.2	MECHANICS.....	5
1.3	ELECTRONICS .....	5
<b>2</b>	<b>QUICK START .....</b>	<b>7</b>
2.1	SETUP .....	7
2.2	COLLECTING DATA .....	7
2.3	DOWNLOADING LOGGED DATA.....	8
2.4	DISCONNECTING.....	8
2.5	SWITCH-ACTIVATED LOGGING.....	8
<b>3</b>	<b>OPERATIONAL DETAILS .....</b>	<b>9</b>
3.1	POWER AND DATA INTERFACE .....	9
3.2	DATA SAMPLING AND LOGGING.....	10
3.3	POWER-UP, SLEEP MODE AND SWITCH CONTROL .....	10
3.4	DEPLOYMENT .....	12
3.5	CLEANING .....	14
<b>4</b>	<b>BATTERY PACK.....</b>	<b>15</b>
4.1	DESCRIPTION.....	15
4.2	CHARGING.....	16
<b>5</b>	<b>FIRMWARE COMMANDS.....</b>	<b>17</b>
5.1	PROTOCOL.....	17
5.2	COMMON COMMANDS.....	17
<b>6</b>	<b>DATA PROCESSING REFERENCE.....</b>	<b>21</b>
6.1	DATA FORMATS .....	21
6.2	.RAW FILE FORMAT .....	21
6.3	.DAT FILE FORMAT.....	22
6.4	CALIBRATION FORMULAS.....	23
6.5	CALIBRATION FILE FORMAT.....	24
<b>7</b>	<b>MAINTENANCE, SERVICE AND SUPPORT .....</b>	<b>27</b>
7.1	MAINTENANCE .....	27
7.2	SERVICE AND SUPPORT.....	27
7.3	EQUIPMENT RETURN .....	27

# 1 INTRODUCTION

The Gamma-2 Abyss is a dual-wavelength transmissometer that operates at ocean depths up to 6,000 meters. It is part of the HOBI Labs Gamma family of spectral transmissometers, which differ in depth ratings and number of wavelengths, but share many features. The Gammas use a folded optical path to keep all electronics in a single pressure housing, and to support multiple path length options up to a full meter (the standard path length is 30 cm). The optical path is completely open to natural water flow so no pumping is required, and the instrument contains no moving parts. The electronics include an intelligent controller with internal logging and large flash memory capacity.

## 1.1 Opto-electronic Design

Figure 1 shows a schematic view of the beam optics. The two wavelengths are generated by separate LEDs. Each LED's current is controlled through software, and modulated at a frequency and phase that allows it to be distinguished from the other. The LEDs' beams pass through diffusers (not shown) and then a beam splitter that directs a portion of each to the output beam, and a portion to a reference photodiode. The reference photodiode and its associated preamplifier electronics are identical to those used to receive the signal after it passes through the water, except that the signal preamplifier has much higher gain. For both the reference and signal measurements, the signals from the two LEDs are separated from each other through digital synchronous detection.

The source and receiver lenses have 100 mm focal lengths, and 1 mm pinholes at their focal planes make the beam divergence and receiver field of view less than  $0.5^\circ$  full-width in water. The beam exits the source window with a 5 mm diameter, and over a 1 m path expands to 14 mm.

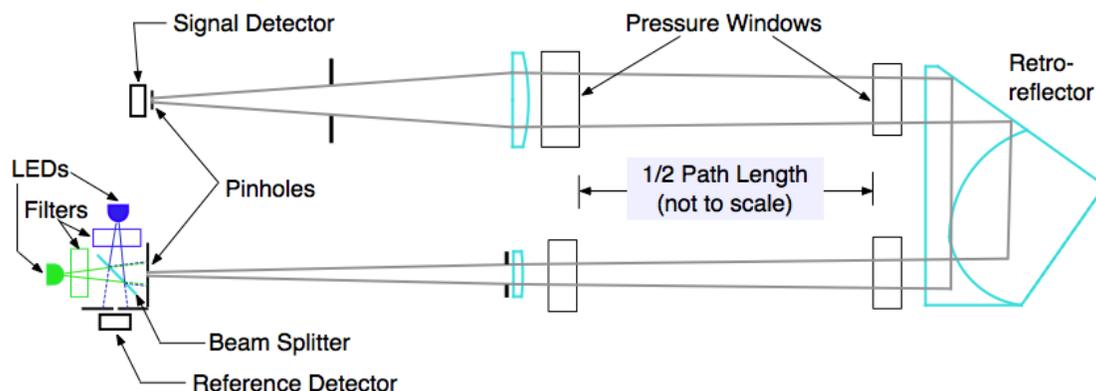


Figure 1. Optical Schematic

Note that ideal attenuation measurements require optics that completely reject all scattered light, regardless of how small the scattering angle. In practice, all transmissometers have finite acceptance angles, receive some scattered light, and therefore report attenuation values smaller than the theoretical value. Gamma-2 Abyss's optics are closer to the ideal than most transmissometers, and therefore it will report somewhat higher attenuation values than those instruments. For the same reason, it is also more sensitive to individual particles and other optical perturbations, so it reveals more of the inherent variability in the properties of natural water.


## 1.2 Mechanics

The Gamma-2 Abyss achieves its extreme depth rating through the use of titanium for the housing, and custom-made sapphire windows. The retroreflector is isolated from ambient pressure in its own housing with sapphire windows. The path length is determined by the standoff that attaches the retroreflector housing to the main housing. The standoff is machined in one piece to give it great strength and rigidity. Figure 2 shows the instrument's outside dimensions.

## 1.3 Electronics

A powerful logging computer, 128 MBytes (more upon request) of nonvolatile flash memory, and integration with a companion battery pack (section 4) allow the Gamma-2 Abyss to operate autonomously for deep profiling or long-term sampling, without real-time electrical connections to the surface.

A high-resolution multi-channel digital-to-analog converter sets the LED currents under the control of the computer, allowing the current to be trimmed for optimum signal levels and balance between the wavelengths. The currents to the LEDs are modulated and synchronized with the analog-to-digital conversions so that the signals from the LEDs can be digitally isolated from each other and processed separately, even though they share detectors and

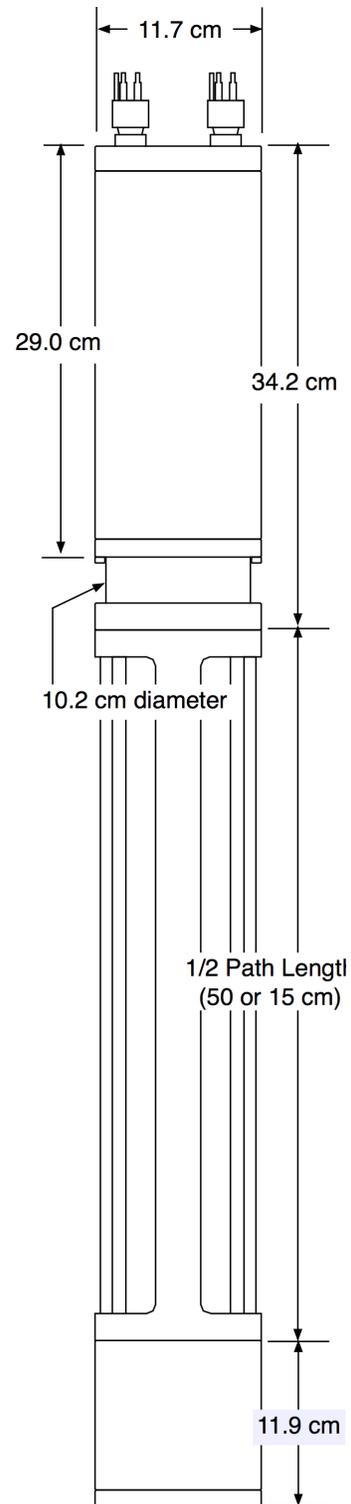

**Figure 2. Outside Dimensions**

preamplifier electronics. All signals are digitized at 1000 samples per second and digital filtering is applied to minimize noise in the final outputs.

## 2 QUICK START

### 2.1 Setup

- Install HydroSoft on your computer. HydroSoft 2.80 or later is required.
- Connect the power/data interface (PDI) cable to the appropriate RS-232 serial port on your computer, or to a USB-Serial adapter.
- Connect the PDI's underwater connector to the Gamma-2 Abyss.
- Run HydroSoft by selecting it on the Windows **Start** menu.
- Click  (or select **Connect** from the **Instrument Menu**).
- Click the **Search** button.
- After a short time HydroSoft should identify the instrument, load its calibration, and close the **Connect** dialog box.
- Select **Set Date/Time...** from the **Gamma-2** menu.
- Click **Set Time** to synchronize the Gamma-2 clock with your computer's clock. This will take a few seconds. Then click **Close** to close the **Date/Time** dialog box.

### 2.2 Collecting Data

- Select **Sampling Options...** from the **Gamma-2** menu.
- For testing purposes, we recommend the following settings:
  - Sample Rate: 2 per second
  - Start Delay: 0 seconds
  - Burst mode: unchecked
- Select **Plot Vs. Time** from the **Graph** menu.
- Click  (or select **Start** from the **Gamma-2** menu). This will cause data to be logged in the flash memory as well as sent to HydroSoft.
- If no data are visible after a short time, click , and the plot will zoom as needed to show all data.
- To adjust the way data are displayed, double-click on the graph, or select **Properties...** from the **Graph** menu.

- Click  (or select **Stop** from the **Gamma-2** menu).
- Close the open data file by selecting **Close** from the **File** menu, or clicking on the  button. If desired, save the file when prompted.

### 2.3 Downloading Logged Data

- Select **Get Data From Instrument** from the **Gamma-2** menu.
- After several seconds a directory of casts in memory will appear. The last cast listed should be the one you just recorded. If that cast is not already highlighted, click on it to highlight it, or click the **Select All** button.
- Click the **Browse...** button to select a directory in which to save downloaded data files. In the resulting dialog box, select a directory and click the **Select** button.
- If it is not already checked, check the **Create calibrated data (.dat) files** option.
- Click the **Download** button to begin transferring data casts.

### 2.4 Disconnecting

- After you are finished communicating with the **Gamma-2**, select **Disconnect** from the **Gamma-2** menu. When asked “Put **Gamma-2** to sleep before disconnecting?” click **Yes**.
- Unplug the PDI cable from the **Gamma-2**.

### 2.5 Switch-Activated Logging

Note that the switch behavior can be modified with firmware commands. For details see section 3.3.

- To start a cast when the **Gamma-2** is asleep and disconnected from the computer, move the switch on the end cap to the start position (in the direction of the arrows on the switch actuator). If the switch is already in the start position, first move it briefly to the opposite position.
- To end a cast, return the switch to the stop position.

## 3 OPERATIONAL DETAILS

### 3.1 Power and Data Interface

#### **Power Supply**

The Gamma-2 requires a power source of 10 to 24 V at 1 to 1.5 W. The exact power consumption depends somewhat on the choice of wavelengths, since different LEDs have different efficiencies. With 470 nm and 532 nm channels, power is typically 1.2 to 1.3 W over the entire supply voltage range.

In sleep mode (see section 3.3), current consumption is less than 100  $\mu$ A at 12V.

#### **Connectors**

Data and power connections are through a SubConn BH4M connector (mating connector: IL4F). A BH2M connector (mate: IL2F) provides a second power input. Figure 3 shows the connector pin assignments. The power inputs on the two connectors are identical and protected against reverse voltage, so it is safe to connect separate power sources to them. Power will be drawn from the source that has the highest voltage.

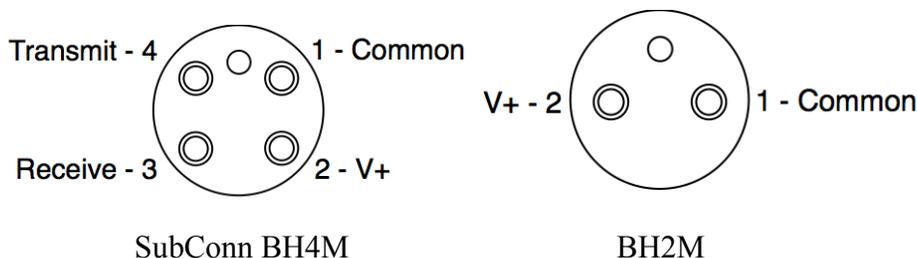


Figure 3. Connector pin diagrams (face view of male connectors)

#### **Power/Data Interface (PDI) Cable**

The optional PDI cable connects to the 4-pin connector and provides 4-way binding posts for power and a DB9F for connection to a computer's RS232 port.

#### **Data interface**

Commands and data are transmitted via RS232 (8 bits, no parity, 1 stop bit). The default baud rate is 57600 baud, with standard rates from 2400 to 115200 supported. In HydroSoft, the baud rate can be changed from the Gamma-2 menu. See section 5 for information about commands, and Section 6 for data formats.

## **3.2 Data Sampling and Logging**

Data are logged in flash memory according to parameters set with HydroSoft's *Sampling Options...* command, or with the LOG firmware command (section 5.2.9). Data casts are started and ended either by software command or by the switch on the end cap (Figure 4). Depending on the LOG and switch settings, a data cast may also be started each time power is applied to the instrument.

During casts, every sample saved in memory is simultaneously transmitted from the serial port.

Individual samples can also be requested with the D firmware command (section 5.2.3). Samples requested this way are not logged.

### ***Averaging and Packet Timing***

Gamma-2 collects and averages data internally at 1000 Hz, and the values transmitted and logged are averages of these rapid samples. In general, a packet is an average of all the values collected since the previous packet, and until the moment the current packet is sent. In other words, if it is logging data once per second, the value reported at time  $X$  is actually the average of the value over the times from  $X-1$  to  $X$  seconds. The time value reported in the packet is  $X$ , the end of the sampling interval. For sample periods up to 8 seconds, all the rapid samples are given equal weight in the average. For longer periods, older samples are given gradually decreasing weight.

At the beginning of a data cast (started with the switch or START command), older samples are flushed from the average so that the first packet reported will include the same number of averaged samples as subsequent packets.

## **3.3 Power-up, Sleep Mode and Switch Control**

Gamma-2 provides various options to control its power mode and its response to the switch on its end cap (Figure 4). For most applications, there is no need to change from the default settings.


### **3.3.1 Default Setup**

The default configuration provides simple control of logging with the switch or with real-time control by HydroSoft. Setting the switch to “start” wakes the instrument, if it is asleep, and starts logging. Setting the switch to “stop” ends the current logging cast and puts the instrument into low-power sleep.

To restore the default settings, send the commands SW,1,1,1 and LOG,,0. To make the LOG settings permanent, also send STORELOG. The SW settings are automatically made permanent without a separate command.

### 3.3.2 Sleep-Wake Details

In sleep mode, all activity is suspended and current consumption falls to about 100  $\mu$ A. Sleep can be initiated either with the external switch or with the SLEEP command (see section 5.2.11). During sleep all settings and state information are maintained in battery-backed memory. If a logging cast is in progress at the time of sleep, it will resume upon the next wake (although if sleep is triggered by the switch, by default the cast will be stopped first).

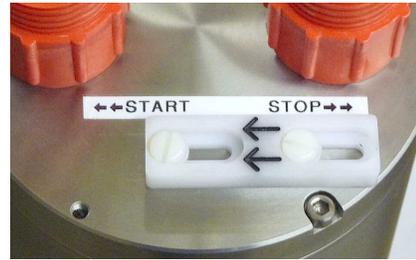

**Figure 4. Switch location. The arrows on the switch actuator point toward the “start” position.**

Sleep is ended by moving the switch to the start position, by any activity on the serial port input, or by the real-time clock if a sleep time was specified with the SLEEP command.

If power is removed while the Gamma-2 is asleep, applying power again will **not** wake it unless the switch is in the start position (see Switch Details below)

### 3.3.3 Switch Details

The behavior of the switch can be altered from the defaults described above with the SW firmware command (also see section 5.2.15). SW can independently set whether the switch starts casts, stops casts, and puts the instrument to sleep. However the switch will *always* wake the Gamma-2 from sleep.

Note that, except at the time power is first applied, the switch only acts when moved from one position to the other, and its position can be overridden by software commands. For example, if the switch is in the stop position you can still start a data cast from HydroSoft or by sending the START command. If you move the switch to the start position while that cast is still in progress, it will have no effect, but moving it back to the stop position will stop the cast (if the *StopsCasts* parameter of the SW command has its default setting).

The one case in which the switch can affect the instrument’s behavior without being moved from one position to the other is when power is first applied. If the switch is in the start position, applying power will always wake the instrument, and it will start a cast if the switch’s *StartsCasts* parameter is set to 1.

### 3.3.4 Power-up Behavior

Figure 5 shows how the Gamma-2 reacts to the application of external power, depending on its various command options. When all options are set to their defaults and the switch is in the “stop” position, the instrument will either boot into the “wait for

commands” state, or, if was put to sleep before power was last removed, it will remain asleep. If the switch is in the “start” position when power is applied, the Gamma-2 will always wake.

If it is not waking from sleep, the Gamma-2 will reboot each time power is applied. The reboot is safe and normal, but resets volatile parameters such as logging settings that have not been stored with the STORELOG command (see section 5.2.14). The real-time clock setting is always retained.

### 3.3.5 Low-Voltage Sleep

To prevent over-discharging external batteries, the Gamma-2 has a programmable low-voltage threshold; if the supply voltage falls below the threshold, the Gamma-2 will go to sleep. You can adjust the threshold, or turn off this function, with the VSLEEP command (see section 3.3.5).

### 3.3.6 Using Power to Control Logging

Note that while it is possible to use the Gamma-2 in a “dumb” mode where you initiate logging simply by applying power to it, you should not routinely end casts by simply removing power. Cutting power while a cast is in progress will definitely result in the loss of up to 30 seconds of buffered data from the cast, and in addition, creates the risk of damaging the flash file directory and losing all logged data completely.

## 3.4 Deployment

### 3.4.1 Mounting

In order for the Gamma-2 to be safely lowered on a cable, it must be secured in a frame or other structure that provides a suitable attachment point. Figure 6 shows the recommended frame, supplied by HOBI Labs. The figure also shows a battery pack attached to the frame, which is required if the lowering cable does not have conductors for supplying power.

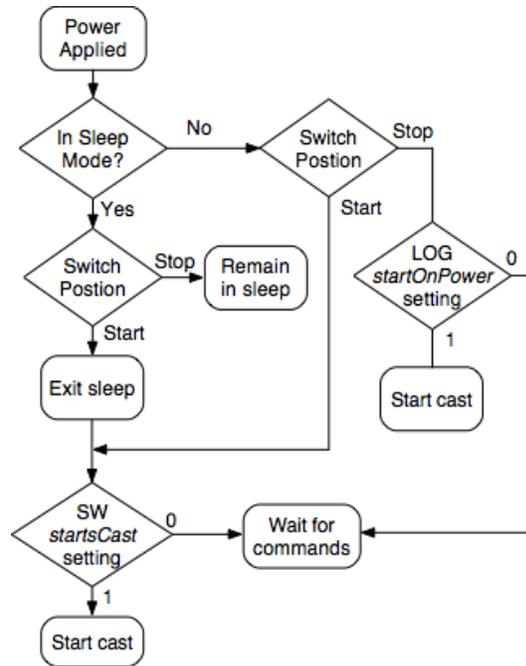


Figure 5. Program flow when power is first applied.

The frame shown here has several features that make it very secure if it is assembled and used properly. It consists of two pairs of clamps that go around the main housing of the Gamma-2. The lower pair has a smaller inside diameter that matches the diameter of the Gamma-2's "waist" and ensures that it cannot possibly slide off the instrument. The upper pair of clamps is rigidly connected to the lower pair with stainless steel rods so the clamp pairs provide mutual support. A forged stainless steel eye is threaded into the top of each rod, and it is to these four eyes that user can secure a lowering harness. If using rope to attach to these eyes, do not use a single loop of rope threaded through all four, because that would be vulnerable to a single-point failure by the rope. Use one loop for each pair of eyes, or some other redundant arrangement.

Figure 6 also shows a clamp pair on the lower end of the instrument. This simply provides convenient protection for the instrument housing during handling, especially when it is laid on its side on deck. It also allows it to stand on end as shown in the figure, but **never** leave the Gamma-2 standing unsecured and unattended.

### 3.4.2 Real-time cable connection

If your deployment cable has the necessary electrical conductors, you can power the Gamma-2 and communicate with it in real time, including controlling it with HydroSoft. However the limitations of typical electrical cables make this impractical for applications that exploit the Gamma-2's extended depth capability.

To provide adequate power to the instrument, the cable's round-trip resistance must be 100 ohms or less, and real-time communication may require even lower resistance. This limits practical cables to at most a few hundred meters length and possibly much less. If you intend to use a real-time cable connection for deployment, always test the cable and instrument beforehand, in the deployment configuration.

### 3.4.3 Battery-powered deployment

The Gamma-2 works seamlessly with a battery pack to provide data logging without any electrical connection to the surface. For battery-power profiling you can



**Figure 6. Gamma-2 Abyss with deployment frame and battery pack (also see section 4).**

leave the battery connected continuously and start and stop logging with the end cap switch. Each switch-activation will create a new numbered cast in flash memory.

### 3.5 Cleaning

Cleanliness of the windows is critical to the Gamma-2's measurements. Because it is sensitive to fractions of a percent change in light intensity, even a tiny speck of contamination can make a significant difference in the readings. Fortunately cleaning is not difficult, and the instrument's sapphire windows are difficult to damage.

The degree of cleaning required depends on the environment in which the instrument is operating, and the deployment routine.

The most basic, but one of the best, cleaning steps is to rinse the windows with generous quantities of distilled or deionized water. You should always finish your cleaning process with a pure-water rinse, and it is highly advisable each time the instrument is withdrawn from the water, so as to prevent substances from drying onto the windows. If the water you are measuring is among the clearest natural waters, this alone is probably sufficient cleaning between casts.

For more thorough cleaning and work in "dirtier" waters, you can use solvents including

- ethanol,
- weak detergent solutions, and
- lens cleaners (such as part number NT54828 at <http://edmundoptics.com>).

For the most thorough possible cleaning, use each of these in the order listed, and follow with a pure-water rinse. In relatively clear water it is adequate to do such a thorough cleaning only once per day. In water rich with growth, sediment, or oily substances, more frequent use of solvents may be necessary, always followed with clean water.

**NOTE:** do not use acetone or other harsh solvents, which can attack the plastic retaining rings and o-rings.

## 4 BATTERY PACK


The Abyss battery pack is an ideal companion for the Gamma-2 Abyss in deep profiling applications. It matches the Gamma-2 Abyss's 6000 meter depth rating and has enough rechargeable capacity to power the Gamma-2 for at least 24 hours of continuous data collection.

### 4.1 Description

The battery pack case is 8.7 cm (3.45") in diameter, 24.3 cm (9.55") long, and is made of hard-anodized 7075 aluminum. The anodizing and a zinc anode on the connector end cap protect against corrosion in water as long as there is no electrical contact between the case and other metals. Plastic mounting brackets such as those supplied by HOBI Labs (see Figure 6) ensure proper insulation.

The pack includes a pressure relief valve to prevent hazardous buildup of gas that could result from (rare) faults in the batteries or charger.

The battery pack has a nominal voltage of 10.8, with the actual voltage ranging up to 13 V when fully charged, depending on temperature. HydroSoft displays the voltage whenever the instrument is connected and logging (but remember that if voltage is being supplied by another source, only the higher voltage is displayed). The battery is roughly 90% discharged when its voltage declines to 10 V. Its total capacity is 3 A-hr, enough to power the Gamma-2 continuously for well over 24 hours.

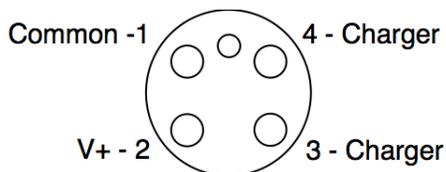

**Figure 8. Battery pack connector wiring. SubConn BH4F connector, face view. Mating connector is IL4M.**



**Figure 7. Abyss Battery Pack**

The output is protected against heavy overload by a positive temperature coefficient thermistor ("resettable fuse") that reduces the output voltage if the current goes above about 3 A. Once tripped, the thermistor will keep the current low (but not zero) until the load is disconnected long enough for the thermistor to cool.

## 4.2 Charging

- **Important:** Charge the pack only with the HOBI-supplied charger.
- Plug the charger into the battery pack **before** connecting the AC input to the charger, and **after** disconnecting the AC input. This prevents shorting the exposed pins on the charger's adapter cable while power is applied.
- Do not charge the battery too frequently. Note that the charger replaces charge 10 to 20 times as fast as the Gamma-2 consumes it.

## 5 FIRMWARE COMMANDS

The Gamma-2 is configured and controlled with commands sent via its serial communication port. When you use HydroSoft, the commands needed for routine operation are generated automatically and you do not need to know the details. However you can type commands manually and view the instrument's replies directly, via HydroSoft's terminal window or a generic terminal program.

### 5.1 Protocol

The Gamma-2 communicates through a standard RS-232 serial connection, with 8 data bits, no parity, one stop bit, and no handshaking. The default baud rate is 57600, but it can be set to other standard rates from 2400 to 115200 through HydroSoft or the BAUD command described in section 5.2.1.

The Gamma-2 waits until it receives a carriage return and/or line feed before responding to a command. When it is not actively logging data, it retransmits each character it receives. Character echo is disabled during logging, so that commands can be sent without affecting the output data stream.

Commands are case-insensitive. They are shown below in upper case for clarity.

Some commands accept arguments, which are separated from the base command, and from each other, by commas or spaces. Arguments may be individually omitted if commas are used as placeholders for them. For example, if a command accepts three arguments, in the form `COMMAND, arg1, arg2, arg3`, you may adjust only *arg3* by entering `COMMAND, , , arg3`.

### 5.2 Common Commands

#### 5.2.1 BAUD *newrate* /S

Reports the current baud rate of the instrument. If a valid *newRate* is specified (2400, 4800, 9600, 19200, 38400, 57600 or 115200), the rate will immediately change, unless the /S switch is included. If /S is included, *newRate* will be saved as the default rate to be used when the instrument is next reset..

Note that, because of the time required to transmit data packets, lower baud rates may limit the maximum sampling rate.

### **5.2.2 CLEARLOG**

Deletes all cast files from flash memory, and resets the cast number to 1. Because this will result in permanent loss of data, the firmware asks for confirmation before proceeding. Deleting large and numerous files can take up to several minutes.

### **5.2.3 D**

Sends a packet with the most recently collected data (averaged over the time since the previous packet was sent), in the format specified by the DATAFORMAT command.

### **5.2.4 DATAFORMAT *formatNumber***

Sets the format of packets that will be sent during logging, cast downloads, and in response to the D command. As of firmware 1.00, *formatNumber* can either be 0 to indicate a “brief” packet or 1 to indicate a “full” packet (see section 6.1 for details).

### **5.2.5 DELCAST *castNumber***

Deletes the given cast (does not prompt for confirmation). Does not affect the current cast counter.

### **5.2.6 DIR /R *startCast endCast***

Displays a directory of the casts in flash memory, including their start times, durations and number of samples. If /R is specified, all casts will be listed in reverse order. If *startCast* and *endCast* are specified, only casts within that range of numbers will be listed. Do not use /R together with *startCast* and *endCast*; specify *endCast* less than *startCast* to display a range of casts in reverse order.

### **5.2.7 DOWNLOAD *castNumber***

Sends the data from the given *castNumber*, if any, in the current data format.

### **5.2.8 ID**

Displays identifying information about the instrument and its configuration. This is automatically used by HydroSoft when it connects to an instrument.

### **5.2.9 LOG *period delay startOnPower***

Displays the current logging settings, and changes them if valid values are provided. *Period* is the number of seconds between samples (which has a minimum value of 0.1). *Delay* is the number of seconds the first sample will be delayed after logging is initiated (not implemented in firmware version 1.00). If *StartOnPower* is set to 1, a cast will start automatically when power is first applied to the instrument. Note that even if *startOnPower* is zero, a cast may be started upon start/stop switch, depending

on its setup. See section 3.2 for switch details. Changes you make to the LOG parameters will not be retained when power is removed from the instrument, unless you use the STORELOG command to make them permanent.

#### **5.2.10 SENDCAL**

Sends the entire contents of the instrument's internal cal file, if any, through the serial port.

#### **5.2.11 SLEEP *seconds***

Sends the instrument into its low-power sleep mode. Logging and any other activities are stopped during sleep. If *seconds* is specified, the sleep will end after that time interval has passed. Otherwise it will sleep until woken by the switch or serial port activity.

#### **5.2.12 START /N**

Starts logging with the current log settings, using the next sequential cast number. If the /N switch is included, data will not be logged in memory, only transmitted through the serial port, and the cast number will remain unchanged.

#### **5.2.13 STOP**

Stops the current logging cast, if any.

#### **5.2.14 STORELOG**

Stores the parameters of the LOG command in nonvolatile memory

#### **5.2.15 SW *StartsCast StopsCast Sleeps***

Displays the current switch parameters, and changes them if valid values are supplied. Nonzero *StartsCast* means that moving the switch from off to on will start logging (as if the START command were issued). Nonzero *StopsCast* means moving the switch from on to off will stop logging (as with STOP). Nonzero *Sleeps* means moving the switch from off to on will put the instrument into an indefinite sleep (as with the SLEEP command). The default setting is 1 for all three parameters. For more about the switch, see section 3.2.

#### **5.2.16 VER**

Displays the installed firmware version.

### **5.2.17 VIN**

Displays the input supply voltage to the instrument. If voltages are applied to both connectors, only the higher one is reported. The voltage reported will be slightly lower than the actual applied to the connector, because it is measured after the voltage passes through a protection diode. This is the same voltage that is reported in full data packets.

### **5.2.18 VSLEEP *voltage***

Display the supply voltage sleep threshold. If the input supply voltage falls below this value, the instrument will go into indefinite low-power sleep. This is intended to protect external batteries from overdischarge. To turn off this function, set the voltage to 0. For more about sleep functions, see section 3.3.

## 6 DATA PROCESSING REFERENCE

### 6.1 Data Formats

Data are transmitted in raw comma-separated decimal ASCII, with the following parameters:

*time,signal1,signal2,reference1,reference2,pressure,temp1,temp2,temp3,Vin,bgnd,smin,smax,rmin,rmax,N*

- **Time** is in seconds since midnight Jan 1, 1970, sometimes called “UNIX time.”
- **Signal1** and **signal2** are raw digital counts measured by the signal detector from the two optical channels.
- **Reference1** and **Reference2** are the signals measured by the reference detector before the beam enters the water.
- **Pressure** is the raw digital reading of the pressure transducer.
- **Temp1**, **temp2** and **temp3** are temperatures measured at three different places within the instrument (not water temperature). The reported values are 100 times the actual temperature in C.
- **Vin** is the supply voltage input, in mV.
- **Bgnd**, **Smin**, **Smax**, **Rmin** and **Rmax** are state-of-health indicators showing the signal levels at different parts of the circuitry.
- **N** is the number of samples averaged to produce the packet. Normally this will be 1000 times the sample period in seconds.

If you use the DATAFORMAT command to set the data format to 0, a briefer version of the above is produced, which is identical except that it ends after temp3. The brief packets are sufficient for producing calibrated data, so you can speed downloads of large data files by selecting this format. The default DATAFORMAT setting is 1, producing the full packet shown above.

### 6.2 .RAW File Format

HydroSoft’s .raw files start with a header similar to the following, but the exact contents or sequence of the header are not guaranteed to remain the same in all versions of HydroSoft.

```
[Header]
HydroSoftVersion=2.80
CreationDate=04/10/10 14:49:43
FileType=raw
DeviceType=Gamma-2
DataSource=G2100100
CalSource=Gamma-2
Serial=G2100100
Config=100
[EndHeader]
```

Following the header is a complete record of every byte received from the instrument. This will include data packets but may include other informational messages. In the following example, the “START” command sent by HydroSoft is echoed by the instrument, which also announces the start and end of the cast:

```
START
Starting cast 6 in 2 seconds.
1274885398.44,1,2,0,0,5588,-1938,-2278,2325,1117,12511,-23402,-
22826,-22015,-21489,504
1274885398.94,1,3,-1,1,5601,-1933,-2273,2325,1117,12511,-23404,-
22827,-22014,-21483,500
1274885401.94,0,-1,-1,0,5592,-1936,-2277,2326,1117,12513,-23399,-
22794,-22012,-21487,500
Stopped cast 6.
```

### 6.3 .DAT File Format

HydroSoft .dat files start with a header similar to the following. As with the raw header, the exact contents and format are not guaranteed, but the [Channels] and [ColumnHeadings] heading will remain the same. The [Channels] heading announces the names of the primary optical data from the instrument, as displayed in the legends of HydroSoft’s graphs, listed one per line. The line following [ColumnHeadings] is a comma-separated list of the names of the parameters listed in the data lines.

```
[Header]
HydroSoftVersion=2.81
CreationDate=04/10/10 14:49:43
FileType=dat
DeviceType=Gamma-2
DataSource=Gamma-2
CalSource=Gamma-2
Serial=G2100100
```

```

Config=100
[Channels]
"c470"
"c532"
[ColumnHeadings]
Time,Depth,c470,c532,IntT

```

```
[Data]
```

After the [Data] heading, lines of comma-separated data will appear, containing the parameters listed under [ColumnHeadings], for example:

```
40324.6180722222,2.38790,0.1032,0.2018,22.83
```

- **Time** is expressed as the number days since the start of January 1, 1900, which is the native format of Microsoft Excel. The time of day, including fractions of seconds, is included in the fractional part of this number.
- **Depth** is in meters.
- **c1** and **c2** are the calibrated attenuation coefficients at the two wavelengths of the Gamma-2.
- **IntT** is the internal temperature of the Gamma-2, in C.

## 6.4 Calibration Formulas

The following formulas support compensation for several possible effects related to temperature and pressure. Depending on the characteristics of the individual instrument, some portions of the processing may be “turned off” by setting the applicable coefficients to zero in the factory calibration file.

### 6.4.1 Depth

Offset-adjusted pressure reading:

$$P(T) = P - P_0 - p(T) + p(T_{P0})$$

where  $P$  is the raw pressure reading direct from the instrument (in counts) and the temperature-related offset is

$$p(T) = k_{p1}T + k_{p2}T^2$$

Calibrated depth in meters of sea water:

$$D = k_{D1}P(T) + k_{D2}(P(T))^2$$

### 6.4.2 Attenuation

Transmission adjusted for temperature and pressure:

$$\tau(T, P) = \left( \frac{S - S_0}{R - R_0} \right) / a_T(T) a_P(P)$$

where

T is the instrument's internal temperature in C,

P is the temperature-adjusted pressure reading  $P(T)$  in digital counts,

$$a_T(T) = \sum_{n=0}^5 k_{Tn} T,$$

and

$$a_P(P) = \begin{cases} 1, & P < P_1 \\ 1 + k_{\tau PX} \frac{(P - P_1)}{P_2 - P_1}, & P_1 \leq P \leq P_2 \\ (1 + k_{\tau PX}) \sum_{n=0}^5 k_{\tau Pn} P, & P > P_2 \end{cases}$$

where  $P_1$  and  $P_2$  are values of  $P(T)$ .

The fully calibrated attenuation coefficient is given by

$$c = \ln \left( \frac{\tau(T_{PW}, 0)}{\tau(T, P)} \right) \frac{1}{L}$$

where

$L$  is the path length in meters

$\tau(T_{PW}, 0)$  is the transmission measured in pure water, at temperature  $T_{PW}$  and atmospheric pressure.

### 6.5 Calibration File Format

The calibration file, whether stored in the instrument or on the computer running HydroSoft, is a plain text file in the following format. The sequence of parameters is not guaranteed, and it is possible additional parameters could be added in future versions. Therefore any software that reads the file should identify parameters by their labels, not

their sequence. Table 1 shows the plain-text names used in the calibration file for the terms in the formulas above.

Any line in the calibration file may also include a comment after the parameter value, marked with a double-slash (//).

[General]

```
DeviceType=Gamma-2    // Saved by HydroSoft 2.80
Serial=C2100000
Label=HOBI
Config=100
MaxDepth=6000
CalTime= 1270080000 (04/01/10)
```

[Depth]

```
kp1=10.215
kp2=-0.1624
P0=1439
TP0=26.68
kD1=0.3619
kD2=0
```

[Attenuation 1]

```
Name=c470
Lambda=470
DeltaLambda=10
L=1.005
S0=-3
R0=0
kT0=0.9926
kT1=0.0031851
kT2=-0.0002447
kT3=4.8421e-6
kT4=0
kT5=0
P1=50
P2=103
kTauPX=0.1188
kTauP0=.66
kTauP1=1.0771e-5
kTauP2=-7.6957e-10
kTauP3=1.1615e-14
kTauP4=0
kTauP5=0
Tau0=1.00167
TPW=20.77
```

[Attenuation 2]

```
<Same format as Attenuation 1>
```

[End]

**Table 1. Parameters as Named in the Calibration File**

<b>Term</b>	<b>Name in Cal File</b>	<b>Description</b>
$k_{p1}, k_{p2}$	kp1, kp2	Coefficients for temperature-related change in pressure offset
$P_0$	P0	Raw pressure reading at atmospheric pressure and temperature $T_{p0}$ .
$T_{p0}$	TP0	Temperature at which $P_0$ was measured.
$k_{D1}, k_{D2}$	kD1, kD2	Coefficients for conversion of digital counts to meters of sea water.
$L$	L	Path length in meters
$S_0$	S0	Raw signal with path blocked
$R_0$	R0	Raw reference with LED disabled
$k_{T0} \dots k_{T5}$	kT0...kT5	Transmission temperature coefficients
$P_1, P_2$	P1, P2	Pressure thresholds, in temperature adjusted counts, for $a_P(P)$
$k_{\tau PX}$	kTauPX	Coefficient of linear pressure correction $a_P(P)$
$k_{\tau P0} \dots$ $k_{\tau P5}$	kTauP0...kTouP 5	Coefficients for polynomial pressure correction in $a_P(P)$
$\tau_0$	Tau0	Temperature-compensated transmission in pure water, at atmospheric pressure.
$T_{PW}$	TPW	Temperature at which $\tau_0$ measured.

## **7 MAINTENANCE, SERVICE AND SUPPORT**

### **7.1 Maintenance**

For best accuracy and to provide the longest life for your Gamma-2, we recommend you return it annually to HOBI Labs for maintenance and calibration.

Cleaning is the most important maintenance for users to perform. Cleaning of the windows, which affects measurement accuracy, is described in section 3.5. Overall cleaning is also important for protecting instruments against corrosion. Rinse the entire instrument thoroughly with fresh water after each use in salt water. Be sure the Gamma-2 is dry before packing it for extended storage.

### **7.2 Service and Support**

For help and advice regarding the operation of your Gamma-2, contact HOBI Labs at [support@hobilabs.com](mailto:support@hobilabs.com), or call our office at (520) 299-2589. You can also download manuals, software and other technical information from our website, [www.hobilabs.com](http://www.hobilabs.com).

### **7.3 Equipment Return**

To return your Gamma-2 to HOBI Labs for any reason, first contact us for a returned material authorization (RMA) number. You can fill out an RMA request form at our web site, or email us directly at [support@hobilabs.com](mailto:support@hobilabs.com). The RMA is very important for tracking your equipment. Please be prepared to provide us with the instrument serial number, the address to which it should be returned after service, and any special requests. Upon issuing the RMA we will send you shipping instructions.

After we receive the equipment, we will inspect and test it, and send a description and cost estimate for any service required. We will not proceed with service until we receive your authorization for the cost.

Always ship your instrument in its original container, which is designed to protect it against rough handling.